

# Branson™ Sonifier Accessories



## Flow-Thru Horn

### Overview

Designed for pharmaceutical research, the flow-thru horn has two inputs or orifices at the non-vibrating nodal point of the horn. A premixed substance is fed through one of the inputs while the horn is ultrasonically activated. Because two inputs are available, two different types of material can be treated simultaneously. The processed solution exits at the tip of the horn. The horn may be used as a continuous-flow device to collect the solution in one large vessel.

### Typical Applications

- Mixing
- Emulsification
- Dispersion
- Formation of mists and fogs

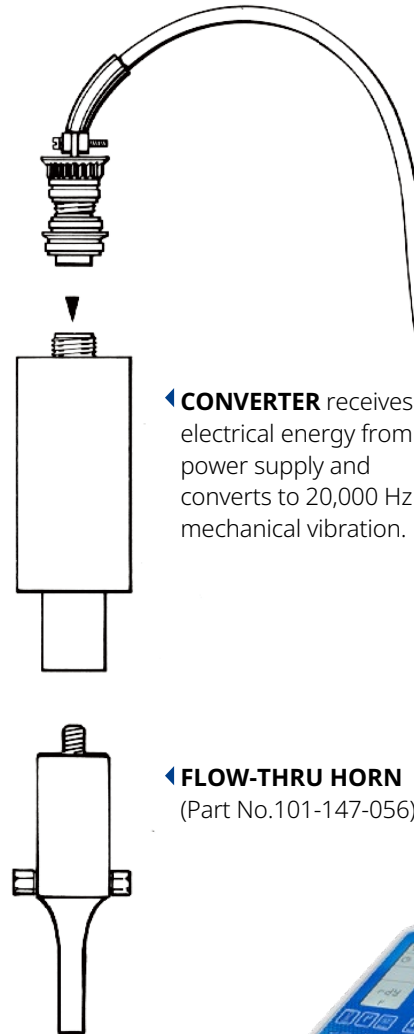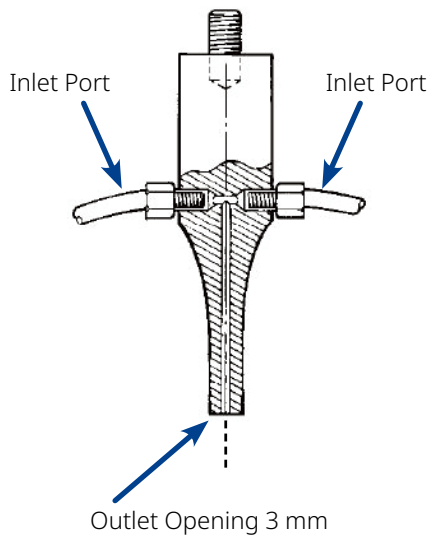

▶ **POWER SUPPLY** converts standard 60 Hz line voltage to high frequency electrical power of 20,000 Hz.

Flow-thru horn compatible with SFX250 and SFX550.



Model Shown SFX550

**BRANSON™** For more information:  
[www.Emerson.com/Branson](http://www.Emerson.com/Branson)

The Emerson logo is a trademark and service mark of Emerson Electric Co. Brand logotype are registered trademarks of one of the Emerson family of companies. All other marks are the property of their respective owners. © 2024 Emerson Electric Co. All rights reserved.

